

PAISEEC

User Manual



Paiseec Travel Mobility Scooter

Simplify Your Travels

www.paiseec.com
support@paiseec.com

CONTENTS

1.Product Introduction	02
1.1 Package Contents	02
1.2 Product Part Names	03
2.Product Functions and Usage	04
2.1 Unfolding the Mobility Scooter	04
2.2 Folding the Mobility Scooter	08
2.3 Carrying the Mobility Scooter	11
2.4 Battery Removal and Installation	12
2.5 Operating the Mobility Scooter	13
2.6 Warning	19
2.7 Indicator Descriptions	22
2.8 Pairing and Removing the NFC Key Card	25
2.9 Voice Prompt Volume Adjust	27
3.Battery and Charging Precautions	29
3.1 Charging Method	29
3.2 Battery Maintenance and Charging Precautions	31
4.Product Daily Care Guidelines	33
5.Product Specifications	34
6.Common Malfunctions and Troubleshooting	36
7.Warranty Terms	37
8.Warnings and Precautions	39

Introduction

Thank you for choosing the Paiseec Travel Mobility Scooter.

The Paiseec Travel Mobility Scooter is an outdoor mobility scooter designed specifically and primarily intended for the elderly and mobility-impaired individuals. We care about your safety, comfort, and ease of use. Please read this manual carefully prior to product use to ensure proper and safe operation.

Please keep this manual stored properly for reference.

The text, illustrations, etc., in this manual may differ slightly from the actual product due to quality improvement or design changes. Our company reserves the right to make product modifications. We appreciate your understanding.

If there are any differences between this manual and the actual product, please refer to the actual product.

For the latest product information, please visit our website (www.paiseec.com).

Please contact us via Email (support@paiseec.com) if you have more questions about our products or services.

Friendly Reminder:


Please read the user manual carefully and familiarize yourself with product performance and operation prior to initial use. Please ensure this manual is properly stored.

1.Product Introduction



1.1 Package Contents

The product package contains the following items:


- Paiseec Travel Mobility Scooter * 1
- Mobility Scooter Key * 2
- Charger * 1
- Quick Installation and User Guide * 1
- NFC Card Key * 2
- User Manual * 1 (this manual)
- Certificate of Conformity * 1

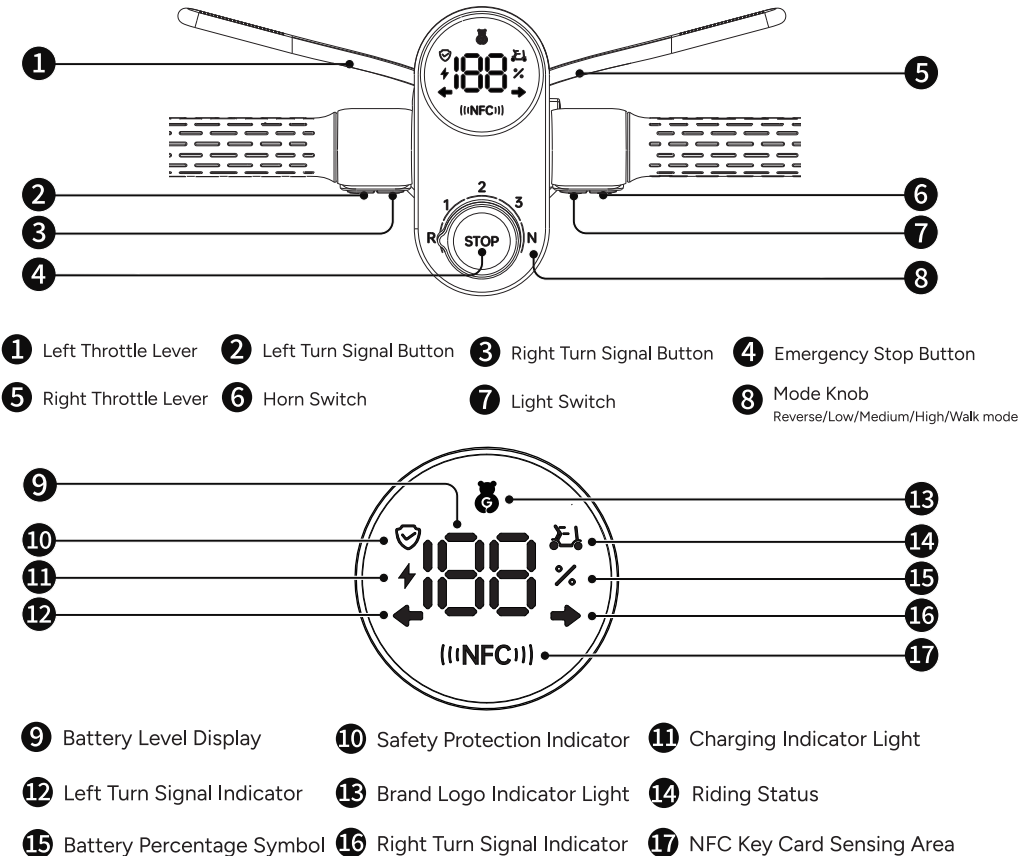

1.2 Product Part Names

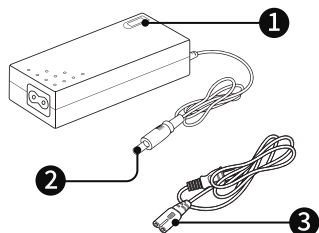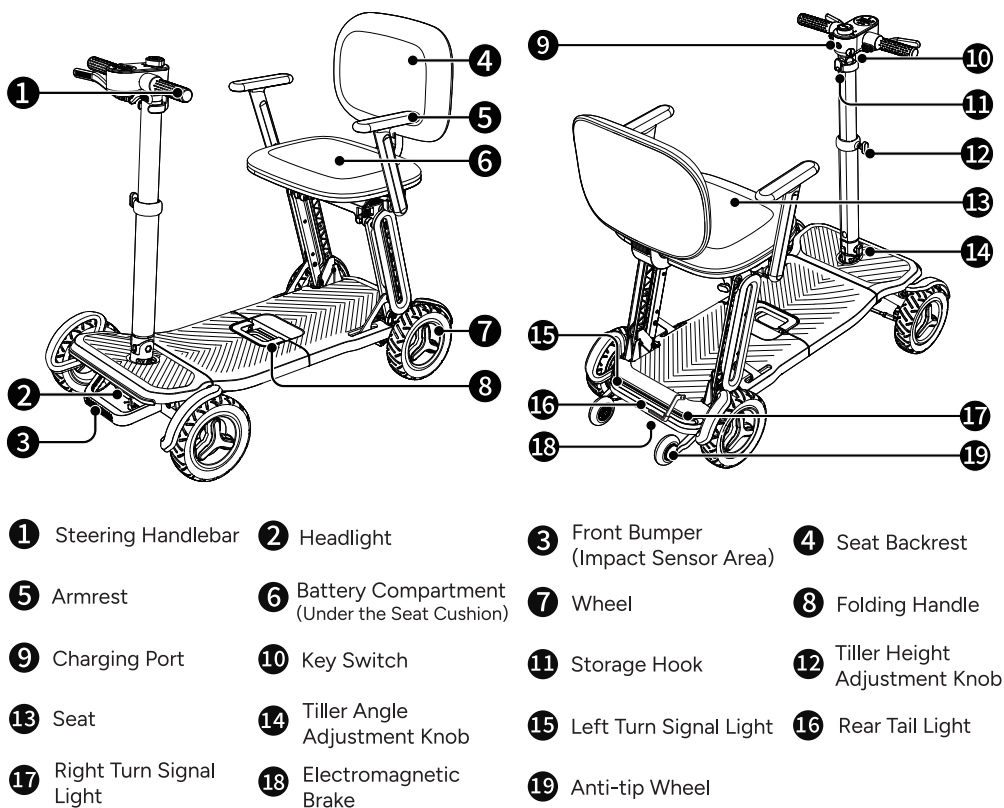

- 1 Charger Indicator
4 Battery Charging Port

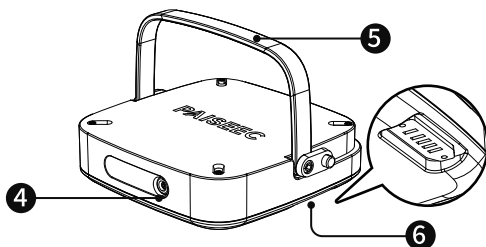

- 2 Charger Connector 3 Power Cord
5 Battery Handle (Battery Release Mechanism) 6 Battery Connector


2.Product Functions and Usage

2.1 Unfolding the Mobility Scooter

- When the scooter is without the seat installed:


1.Unlock the frame:

Push the "Vehicle Unfold/Fold Lock" lever upward from the locked position to the unfold unlock position as shown. Hold the tiller (front steering column) with your left hand and the seat crossbar with your right hand, then pull both ends apart until the scooter is fully unfolded.

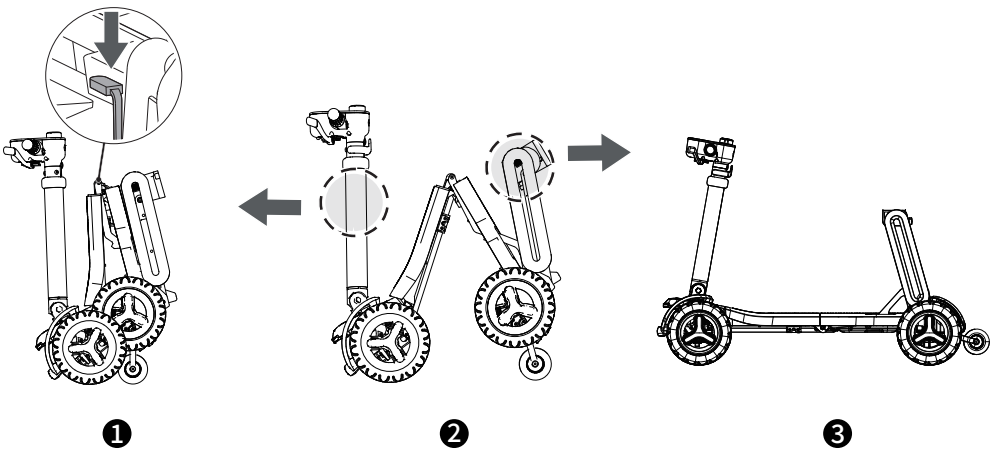


2.Adjust the tiller height:

Rotate the tiller height adjustment knob counterclockwise to loosen it, adjust the tiller to the desired height, then tighten the knob to lock it securely.

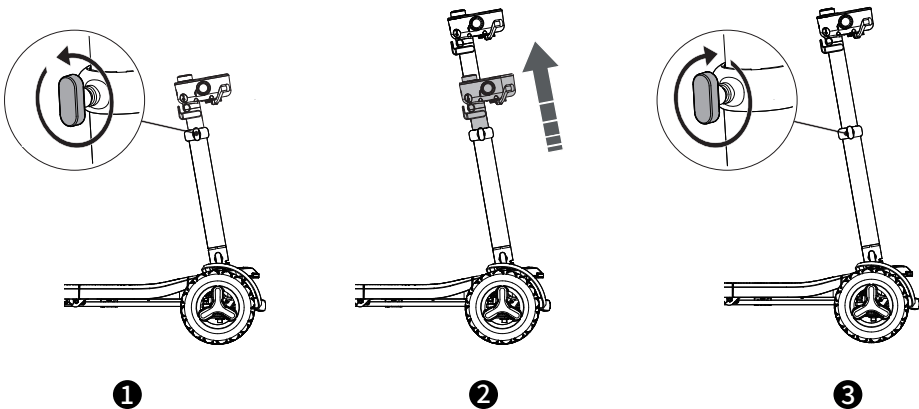


3.Adjust the tiller angle:

Rotate the tiller angle adjustment knob counterclockwise to loosen it, adjust the tiller to a comfortable angle, then tighten the knob to lock it securely.

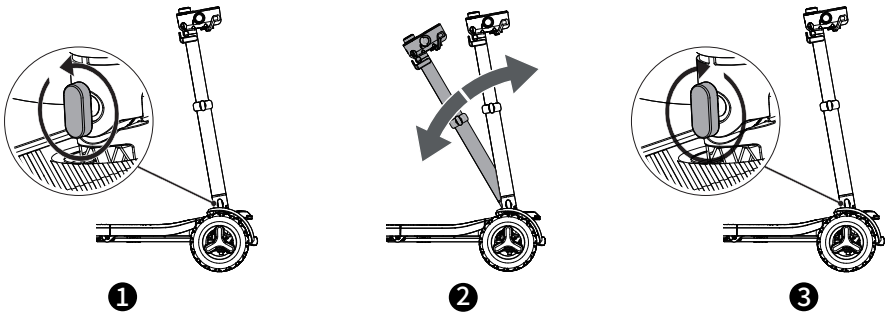


4.Install the seat:

Take out the seat and align the seat post with the mounting hole on the seat crossbar. Insert the post fully, then tighten the seat locking knobs on both sides of the seat. Finally, lift up the seat backrest and armrests into position.

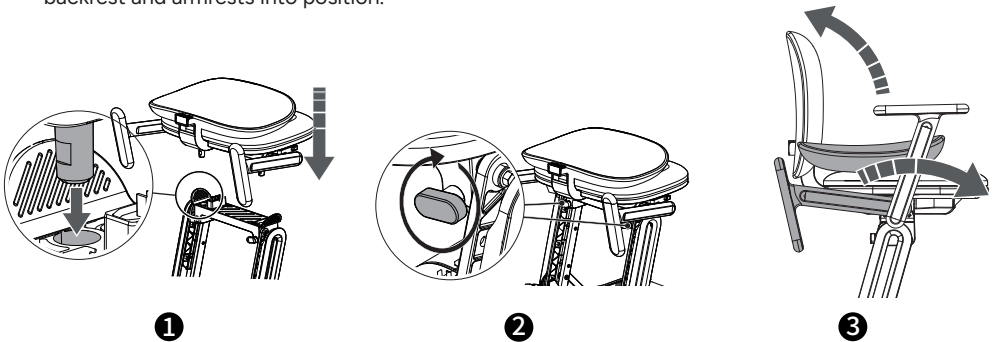

5.Install the storage basket:


Take out the storage basket and open the Velcro straps on both sides. Attach the straps around the slanted support bars on both sides as shown, and press firmly to secure the basket to the chassis. The storage basket installation is now complete.

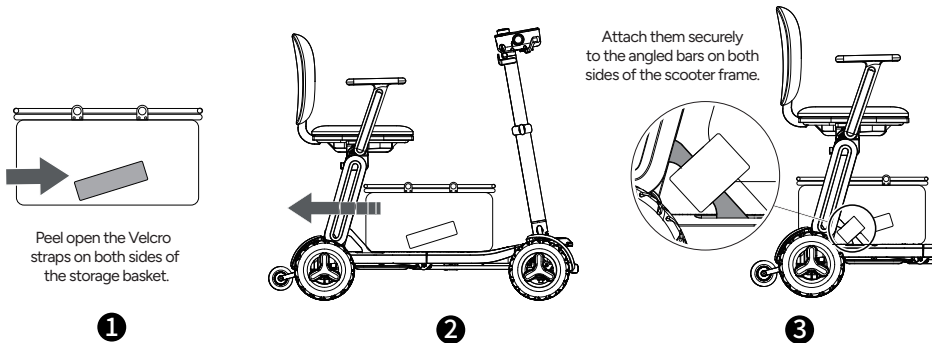

- When the scooter is with the seat installed:


1.Unlock the frame:

Push the "Vehicle Unfold/Fold Lock" lever upward from the locked position to the unfold unlock position as shown. Hold the tiller (front steering column) with your left hand and the seat crossbar with your right hand, then pull both ends apart until the scooter is fully unfolded.

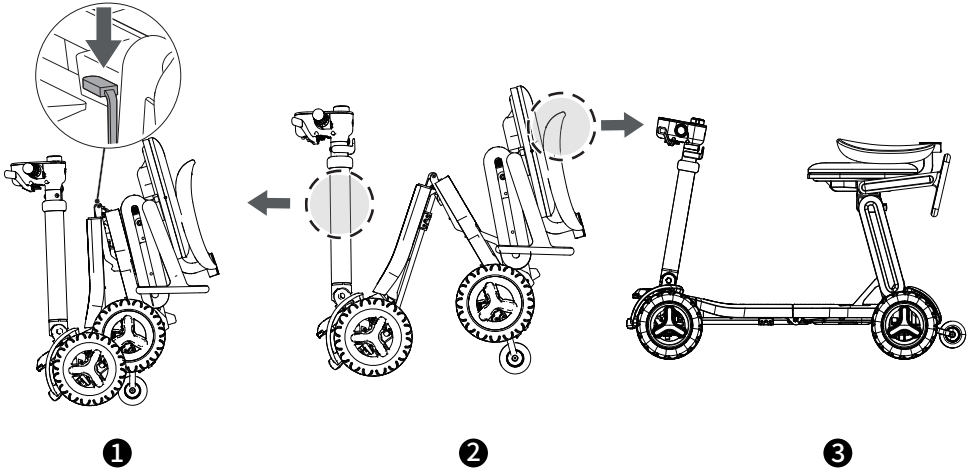


2.Adjust the tiller height:

Rotate the tiller height adjustment knob counterclockwise to loosen it, adjust the tiller to the desired height, then tighten the knob to lock it securely.

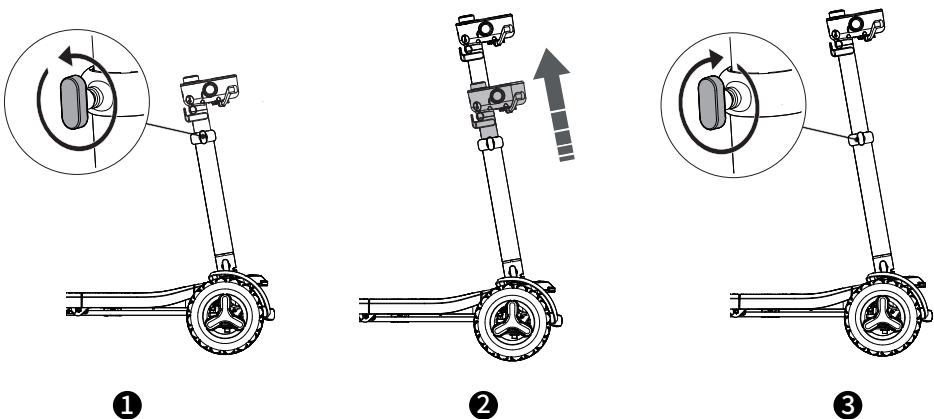


3.Adjust the tiller angle:

Rotate the tiller angle adjustment knob counterclockwise to loosen it, adjust the tiller to a comfortable angle, then tighten the knob to lock it securely.

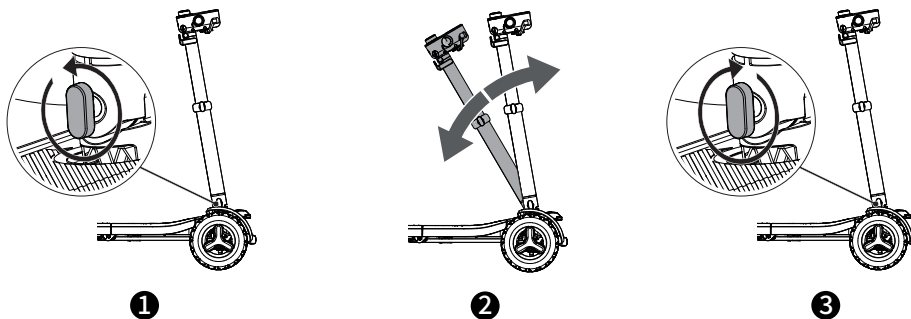


4.Position the seat:

Lift up the seat backrest, then raise both armrests into position.

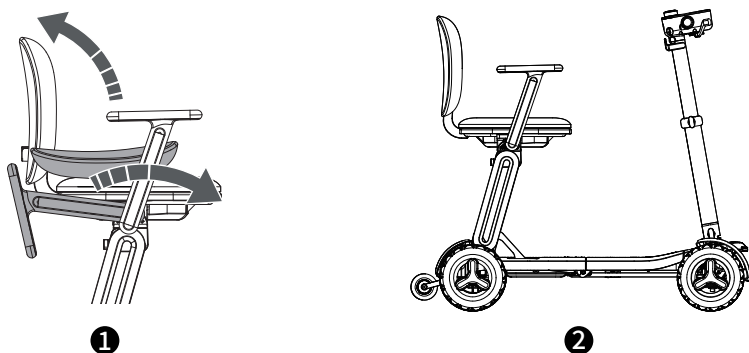

5.Install the storage basket:

Take out the storage basket and open the Velcro straps on both sides. Attach the straps around the slanted support bars on both sides as shown, and press firmly to secure the basket to the chassis. The storage basket installation is now complete.

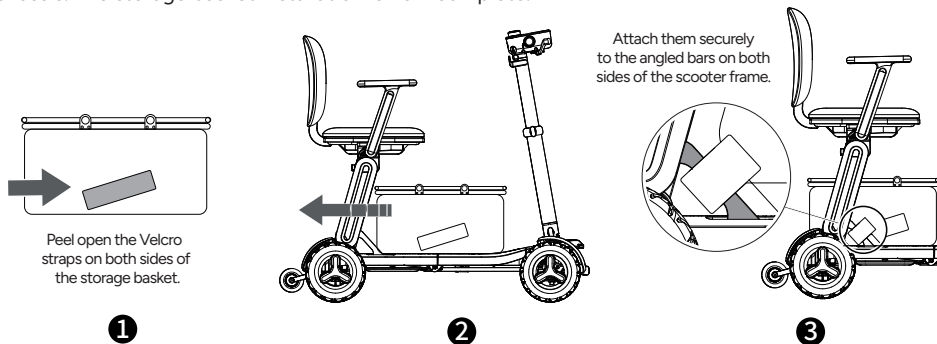


2.2 Folding the Mobility Scooter

- Folding the scooter in the seat-removed state (remove the seat before folding):


1.Remove the storage basket:

Open the Velcro straps that attach the basket to the slanted support bars on both sides, then remove the basket from the scooter chassis.

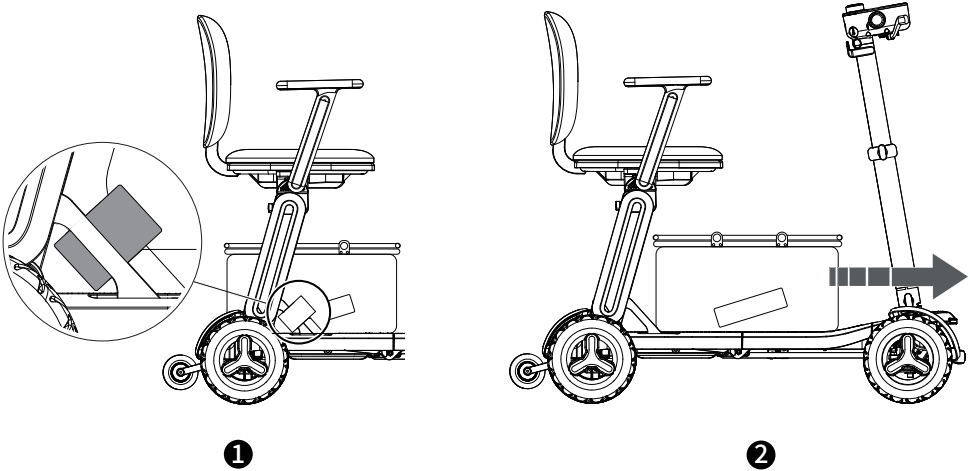


2.Remove the seat:

Rotate the seat locking knobs on both sides under the seat counterclockwise to loosen them. Once loosened, pull the seat post out from the mounting hole on the seat crossbar and remove the seat.

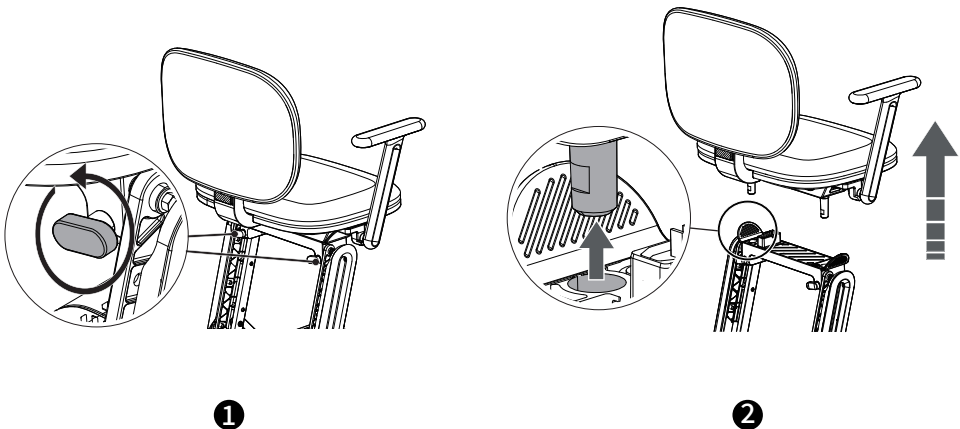

3. Lower the tiller:

Rotate the tiller height adjustment knob counterclockwise to loosen it, then adjust the tiller to the lowest position.

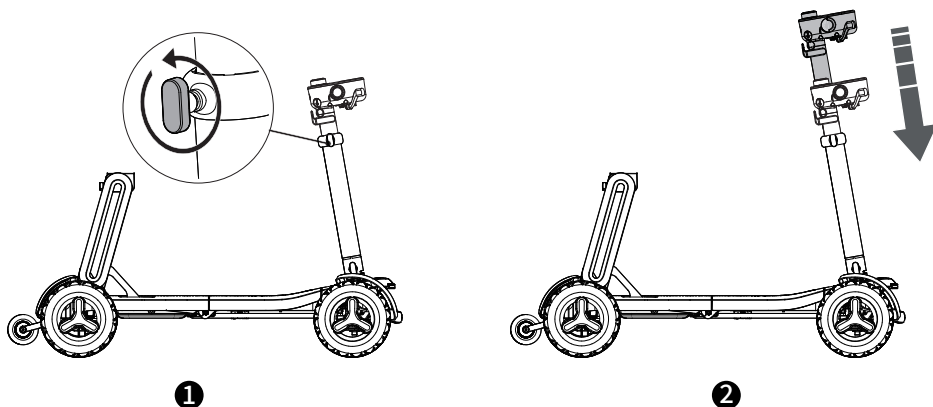

4. Fold the scooter:

Pull up the folding handle to unlock the folding mechanism, then continue lifting the handle to fold the seat, tiller, and chassis inward toward the center. Finally, press the "Vehicle Unfold/Fold Lock" lever downward to the locked position to complete the folding.

